

PANTHERA U3



THE PANTHERA U3 is an active wheelchair designed to let you live life as actively as possible. It is characterized by the attention that we have put into all of its details.

Panthera U3 has a tapered front and a steeper footrest angle that gives a more upright and active driving position compared to the S3.

Its low weight combined with the stable fixed frame and casters with a unique profile on the tread surface make this wheelchair very easy to manoeuvre.

The chassis' seating angle, along with adjustable backrest and seat upholstery and the backrest tilt function combine to provide extremely good seating comfort. The fabric used in the backrest upholstery has also been carefully selected for ergonomic and comfortable sitting.


The chair's balance point can be adjusted to enable you to find a setting that offers a feeling of security. Accessories such as armrests and side guards have a unique design and you can adjust the height and fore-aft position of the armrests yourself without using tools.

STANDARD FEATURES

- Tapered front and steeper footrest angle
- More upright and active driving position
- Chrome molybdenum chassis
- Titanium footbar
- 24" Standard rear wheels
- Balancing, ungraded adjustable
- Carbon fibre rear axle
- unique profiled front castor (3.5" dia) for improved manoeuvrability
- Light weight frame: starting at 9.85 lb (4.48 kg)

COLOURS

PANTHERA S3 Comes standard in a beautiful black frame. It is also available in a series of Premium Colours which can be ordered at an extra cost.

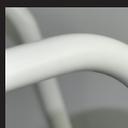

Versatile White



Blue



Mean Green



Panorange



Turtle Turquoise



Brave Bronze

PANTHERA U3

ACCESSORIES


PANTHERA U3 offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately,



Side Guards with a Soft Upper Section which can be used to sit on during transfers or as the bump up to reach higher places such as a cupboard.



Height-Adjustable Armrest Adjusts to 4 heights and 4 horizontal positions.

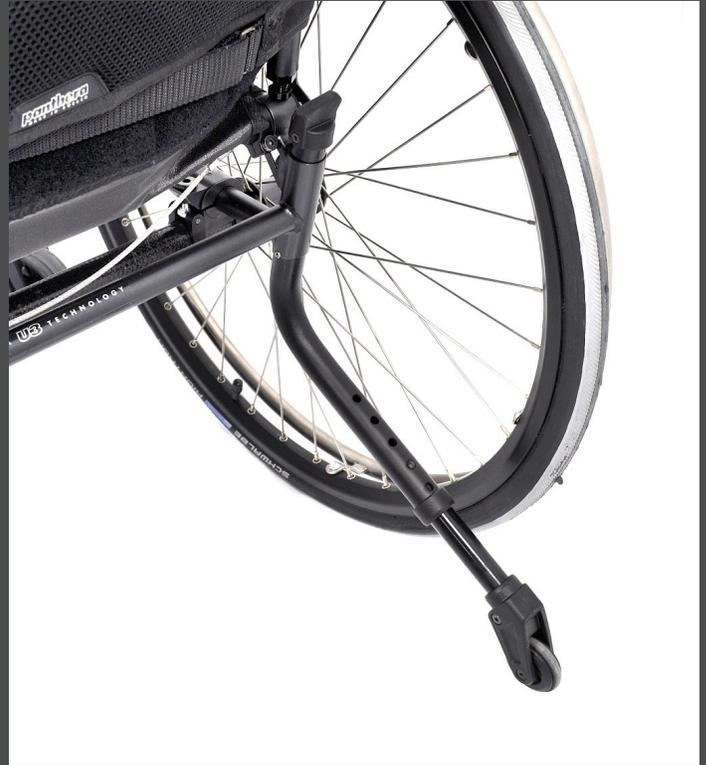

Anti-tip Device for added safety and stability. It can be folded in and out while seated.



Foldable Push Handles



Height-Adjustable Push Handles

PANTHERA U3 SPECIFICATIONS

Seat width	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	18" (45 cm)
Backrest height	8" - 18" (20 - 45 cm)				
Seat height (rear)	17" (43 cm)				
Seat height (front)	18.5" (47cm)				
Seat angle	7°				
Seat depth	13.75" -18.1" (35 - 46cm)				
Total length	30.7" - 35.5" (78 - 90 cm)				
Total height	25" - 33" (64 - 84 cm)				
Total width	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)
Max User weight	220 lb (100 kg)	220 lb (100 kg)	220 lb (100 kg)	220 lb (100 kg)	330 lb (150 kg)
Total weight	17.69lb (8.04 kg)	17.89lb (8.13 kg)	18lb (8.18 kg)	18.08lb (8.22 kg)	18.76lb (8.53 kg)
Transfer weight	9.86lb(4.48 kg)	10.06lb (4.571 kg)	10.25lb (4.66 kg)	10.44lb (4.75 kg)	10.93lb (4.97 kg)



Distributed exclusively
in North America by



Triumph Mobility Inc.

www.triumphmobility.com | info@triumphmobility.com | 1.855.546.0711
R_0617_U3-SalesSheet_EN | Printed in Canada