

# PANTHERA U3 LIGHT



**THE PANTHERA U3 LIGHT** is a wheelchair designed for the experienced active user who can handle a 'tippy' wheelchair without an anti-tip device. After sitting in a wheelchair long term, many people develop problems with neck and shoulder pain. To help reduce and/ or prevent further problems, the U3 light has been designed to make lifting and driving the chair as easy as it can possibly be without, of course reducing the ergonomic features from the U3.

The design is optimized for easy handling in and out of vehicles, and extremely good driving characteristics.

In its standard configuration, the U3 light has a carbon fibre rear axle, Spinerger X rear wheels, X footrest in aluminum and X castor wheels and a seat bag with integrated magnets.

The U3 Light exists in two different configurations. "U3 light" with normal balance, and "U3 light L" with a more "tippy" balance.


On the "U3 Light L" the rear axle has been moved forward an additional 7/8" (22mm). It is possible to adjust the position of the backrest by +/- 0.5" (10 mm) forward or backwards to calibrate the balance.

## STANDARD FEATURES

- Clean & minimalistic design
- Chrome molybdenum chassis
- 'X' footrest
- 'X' 3" (75mm) Castor wheels
- 24" Spinerger X rear wheels
- Carbon fibre rear axle
- One handed brake
- Light weight frame: Starting at 8.26 lb (3.75 kg)

## COLOURS

**PANTHERA U3 LIGHT** comes standard in Pearl White or anodic black. It is also available in a series of premium colours which can be ordered at an extra cost.

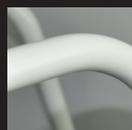


Versatile White



Blue



Mean Green

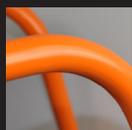

Panorange



Turtle Turquoise



Brave Bronze

# PANTHERA U3 LIGHT

## ACCESSORIES

**PANTHERA U3 LIGHT** offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



**Side Guards with a Soft Upper Section** which can be used to sit on during transfers or as the bump up to reach higher places such as a cupboard.



**Height-Adjustable Armrest**  
Adjusts to 4 heights and 4 horizontal positions.



**Foldable Push Handles**



**Seat Bag**



**Freewheel**

## PANTHERA U3 LIGHT SPECIFICATIONS

|                            |   |                    |                    |                    |                    |
|----------------------------|---|--------------------|--------------------|--------------------|--------------------|
| <b>Seat width</b>          | 13" (33 cm)                                       | 14" (36 cm)        | 15" (39 cm)        | 16.5" (42 cm)      | 18" (45 cm)        |
| <b>Backrest height</b>     | 9" (22 cm), 10" (25 cm), 12" (30 cm), 14" (35 cm) |                    |                    |                    |                    |
| <b>Seat height (rear)</b>  | 17" (43 cm)                                       |                    |                    |                    |                    |
| <b>Seat height (front)</b> | 18.5" (47 cm)                                     |                    |                    |                    |                    |
| <b>Seat angle</b>          | 7°  |                    |                    |                    |                    |
| <b>Seat depth</b>          | 13.8" - 18" (35 cm - 46 cm)                       |                    |                    |                    |                    |
| <b>Total length</b>        | 32.5" ( 82.5 cm)                                  |                    |                    |                    |                    |
| <b>Total height</b>        | 25" - 33" (64 cm - 84 cm)                         |                    |                    |                    |                    |
| <b>Total width</b>         | 21.25" (54 cm)                                    | 22.5" (57 cm)      | 23.6" (60 cm)      | 24.8" (63 cm)      | 26" (66 cm)        |
| <b>Max User weight</b>     | 220 lb (100 kg)                                   | 220 lb (100 kg)    | 220 lb (100 kg)    | 220 lb (100 kg)    | 220 lb (100 kg)    |
| <b>Total Weight</b>        | 13.44 lb (6.1 kg)                                 | 13.53 lb (6.14 kg) | 13.66 lb (6.20 kg) | 13.77 lb (6.25 kg) | 13.88 lb (6.30 kg) |
| <b>Transfer Weight</b>     | 8.26 lb (3.75 kg)                                 | 8.38 lb (3.8 kg)   | 8.48 lb (3.85 kg)  | 8.59 lb (3.90 kg)  | 8.70 lb (3.95 kg)  |



Distributed exclusively  
in North America by



**Triumph Mobility Inc.**

www.triumphmobility.com | info@triumphmobility.com | 1.855.546.0711  
R\_0617\_U3Light-SalesSheet\_EN | Printed in Canada