

# Children's

RIGID WHEELCHAIRS



*panthera*



# TABLE OF CONTENTS

<b>Panthera Bambino 3</b> .....	3
Why choose Panthera Bambino 3 .....	4
Features and Optional Accessories .....	5
Chassis Colour and Specifications .....	6
<b>Panthera Bambino</b> .....	7
Why choose Panthera Bambino 3 .....	8
Features and Optional Accessories .....	9
Chassis Colour and Specifications .....	10
<b>Panthera Micro</b> .....	11
Why choose Panthera Micro .....	12
Features and Optional Accessories .....	13
Chassis Colour and Specifications .....	14



# Why Choose Panthera Bambino 3?

**L**ike Triumph, our partner Panthera is always thinking about ways to innovate. The Bambino 3 children's wheelchair reflects that attitude. Already an easy-to-drive, everyday chair for children with many different needs, Panthera's recent additions offer even more customization suited to children from ages 4 to 12

Here are just some of the benefits built into the Bambino 3:

**Footrest** flips up, moves forward and back, and adjusts vertically for shorter legs, as well as by angle.

**Brake** sits beside the wheel, and is easily accessed by even the most frail children from a lever inside the chair that moves out of the way for transfer.

**Push bar** that's sturdy, quick release and easily adjustable in height and angle. The push bar is removeable to minimize weight in transport and when the user pushes themselves.

**Seat width adaptors** narrow the seat width by 7.5 mm on each side, growing with your child.

**Easily adjustable arm supports, lateral supports and side guards** add customization & comfort.

**Head support** that's cushioned and easily adjustable to the push bar and for multiple head positions for the best comfort and fit.

**Tray table** slides into existing slots for easy use, whether for snacks or browsing on a tablet. An impressive magnet feature means you easily attach it to the wheel of the chair have it to go for the whole day!

**Entire chair easily transported** and fits in the trunk of a car.



## Features


- Lightweight frame: starting at 10.758 lb (4.89 kg)
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- Rear wheels from 20 to 24 inches and 5-inch front caster wheel
- Stable, fixed frame for stability and adjustable balance point

## Optional Accessories

Panthera Bambino 3 offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



**Spoke protector**

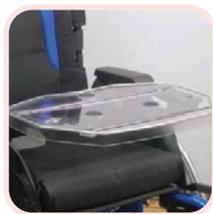

**Tray table**



**Push bar**



**Sideguard**



**Cushion**

# Chassis Colour




Red



Pearl White




Classic Blue

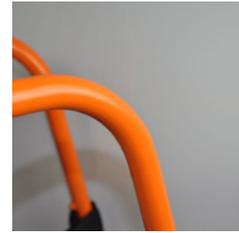

Panorange



Grey Metallic

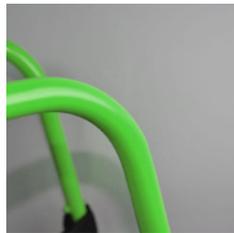

Mean Green



Brave Bronze



Anodic Black

# Specifications

<b>Seat width</b>	9.45" (24 cm)	10.63" (27 cm)	9.45" (24 cm)	10.63" (27 cm)	11.8" (30 cm)	13" (33 cm)
<b>Backrest height</b>	7.87" - 11" (20 - 28 cm) or 10.6" - 13.8" (27 - 35 cm)					
<b>Seat height (rear)</b>	15" (39 cm)					
<b>Seat height (front)</b>	16.1" (41 cm)					
<b>Seat angle</b>	4°					
<b>Seat depth</b>	7" - 9.8" (18 - 25 cm)		9.8" - 11.8" (25 - 30 cm)			
<b>Total length</b>	25.2" (65 cm)		25.6" - 28.1" (65 - 71.5 cm)			
<b>Total height</b>	22.6" - 25.6" (57.5 - 65 cm)					
<b>Total width</b>	21.65" (55 cm)	22.8" (58 cm)	21.65" (55 cm)		24" (61 cm)	25.2" (64 cm)
<b>Max User weight</b>	132 lb (60 kg)					
<b>Total weight</b>	21.98 lb (9.92 kg)				22.48 lb (10.22 kg)	
<b>Transfer weight</b>	12.26 lb (5.57 kg)		12.76 lb (5.80 kg)	13.26 lb (6.03 kg)	12.3 lb (5.6 kg)	



# Why Choose Panthera Bambino ?

**T**he Panthera Bambino offers an easy-to-drive, everyday chair for children with many different needs. The chair offers your child a smooth, stable ride that can raise their confidence about their own mobility and their ability to engage with their busy world.

The chair is ideal for children  
**ages 4 to 12.**

Like the adult versions, the  
Bambino provides a **smooth ride,**  
a **wide range of accessories.**

**Adjustment options** that allow  
it to grow with the user.

The chair is **easy to  
drive, adapt, and lift,**  
and **an abduction  
front allows flexible  
leg positioning,**  
offering **6 cm  
additional width** at the  
knees with the same  
seat depth.



## Features

- Lightweight frame: starting at 10.76 lb (4.89 kg)
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- Rear wheels from 20" to 24" and 5" front caster wheel
- Footrest adjustable in height, angle and depth
- Stable, fixed frame for stability and adjustable balance point
- Dual adjustable antitippers

## Optional Accessories

Panthera Bambino 3 offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



**Spoke protector**



**Push handle**



**Push bar**



**Sideguard**



**Cushion**

# Chassis Colour




Red



Pearl White




Classic Blue

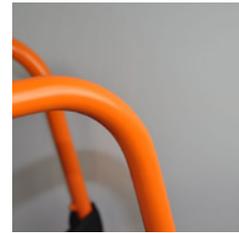

Panorange



Grey Metallic

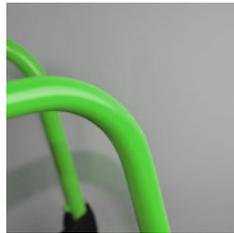

Mean Green



Brave Bronze



Anodic Black

# Specifications

<b>Seat width</b>	9.45" (24 cm)	10.63" (27 cm)	11.8" (30 cm)	13" (33 cm)
<b>Backrest height</b>	6.7 - 9.8" (17 - 25cm) or 10.6 - 13.8" (27-35 cm)			
<b>Seat height (rear)</b>	15" (39 cm)			
<b>Seat height (front)</b>	16.1" (41 cm)			
<b>Seat angle</b>	4°			
<b>Seat depth</b>	7 - 9.8" (18 - 25 cm)		9 - 11.8" (23 - 30 cm)	
<b>Total length</b>	25.6" - 28.1" (65 - 71.5 cm)			
<b>Total height</b>	22.6 - 25.6" (57.5 - 65 cm)			
<b>Total width</b>	20.5" (52 cm)	21.65" (55 cm)	22.8" (58 cm)	24" (61 cm)
<b>Max User weight</b>	132 lb (60 kg)			
<b>Total weight</b>	17.9 lb (8.12 kg)	18.4 lb (8.35 kg)	18.9 lb (8.58 kg)	19.4 lb (8.80 kg)
<b>Transfer weight</b>	10.8 lb (4.89 kg)	11.3 lb (5.14 kg)	11.8 lb (5.37 kg)	12.3 lb (5.6 kg)



# Why Choose Panthera Micro ?

**T**he Panthera Micro is a wheelchair for very young children, designed to meet the needs of those between **1 and 5 years of age** (the youngest user we have met was 8 months old).

**A lightweight** rigid wheelchair for very young children.

Like the adult versions, the Micro provides a **smooth ride** and a wide range of accessories.

The chair is intended for **indoor use, and is extremely easy to drive** for little people who want to be on the move!

The chair is easy to drive, adapt and lift (weight starts at 5.15 lb / 2.342 kg), and offers **16” or 18” rear wheels** and a **3” single front caster wheel** for easy rotation.

**Seat depth adjustment** options and growth pad accessories allow the chair to grow with the user.



## Features

- Lightweight frame: starting at 5.15 lb (2.342 kg)
- Chassis made from chrome molybdenum steel, rear axle from titanium (10° camber)
- 16" or 18" rear wheels and 3" front caster wheel
- Height-adjustable footplate
- No pushrims
- Anti-tipper

## Optional Accessories

Panthera Bambino 3 offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



**Spoke protector**



**Push handle**



**Extended Backrest**



**Sideguard**




**Cushion**


# Chassis Colour




Classic Blue

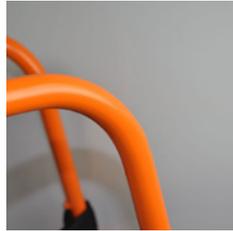

Panorange

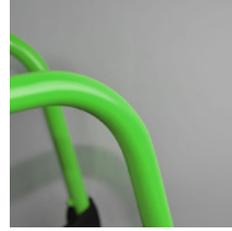

Mean Green

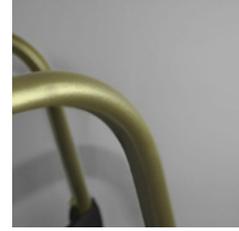

Brave Bronze

# Specifications

<b>Seat width</b>	9.4" (24 cm)	10.6" (27 cm)	9.4" (24 cm) Long	10.6" (27 cm) Long
<b>Backrest height</b>	8.7" (22 cm)			
<b>Seat height (rear)</b>	12.2" (31 cm)			
<b>Seat height (front)</b>	12.1" (31 cm)			
<b>Seat angle</b>	0°			
<b>Seat depth</b>	5.8" - 7.8" (15 - 20 cm)		7.8" - 9.8" (20 - 25 cm)	
<b>Total length</b>	24.4" (62 cm)		28.5" (72.5 cm)	
<b>Total height</b>	20.25" (51.5 cm)		21.25" (54 cm)	
<b>Total width</b>	18" (46 cm)	19.25" (49 cm)	18" (46 cm)	19.5" (49 cm)
<b>Max User weight</b>	66.14 lb (30 kg)			
<b>Total weight</b>	7.76 lb (3.52 kg)	7.86 lb (3.56 kg)	8.34 lb (3.78 kg)	8.44 lb (3.83 kg)
<b>Transfer weight</b>	5.16 lb (2.34 kg)	5.3 lb (2.38 kg)	5.45 lb (2.47 kg)	5.56 lb (2.52 kg)



**PANTHERA EXCLUSIVE DISTRIBUTOR IN NORTH AMERICA**

 **triumphmobilty.com**

 **info@triumphmobility.com**

 **1.855.546.0711**

 **866.225.3059**