

PANTHERA S3 SWING FAMILY



THE PANTHERA S3 SWING FAMILY WHEELCHAIRS have the same great features of the S3 Wheelchairs, but with divided swing away footrests.

The locking mechanism of the footrests is discrete and practical and is integrated inside the frame with a user friendly design.

The S3 Swing has a seat angle of 4 degrees which opens up new possibilities to find a good seating position, even for those who need to use their legs to move about with the wheelchair or for others who need to have a more flat seat angle.

The low weight, combined with stable fixed frames and casters, with a unique profile on the tread surface, make this wheelchair very easy to maneuver. The chassis' seating angle, adjustable backrest, and seat upholstery, with the backrest tilt function, combine to provide extremely good seating comfort.

The chairs' balance point can be adjusted to enable you to find a setting that offers a feeling of security. Accessories such as armrests and side guards have a unique design and you can adjust the height and fore-aft position of the armrests yourself without using tools.


The S3 Swing is available in 2 lengths: S3 Swing and S3 Swing short. The S3 Swing short is 2" (50mm) shorter than the S3 Swing.

STANDARD FEATURES

- Chrome molybdenum chassis
- Folding footrests
- 4" degree Seat angle
- 24" Standard rear wheels
- Carbon fibre rear axle
- Balancing continuously adjustable
- Integrated anti-tip device
- Light weight frame: Starting at 9.15lbs (4.160kg)

COLOURS

PANTHERA S3 SWING comes standard in anodic black. It is also available in a series of premium colours which can be ordered at an extra cost.

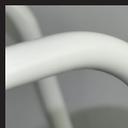

Versatile White



Blue



Mean Green



Panorange



Turtle Turquoise



Brave Bronze

PANTHERA S3 SWING FAMILY

ACCESSORIES

PANTHERA S3 SWING offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Sideguards with a soft upper section which can be used to sit on during transfers or as the bump up to reach higher places such as a cupboard.



Height-adjustable armrest. Adjusts to 4 heights and 4 horizontal positions.



Foldable push handles



Push bar



Seat bag

PANTHERA S3 SWING SPECIFICATIONS

Seat width	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	18" (45 cm)	20" (50 cm)
Backrest height	8" (20 cm), 10" (25 cm), 12" (30 cm), 14" (35 cm), 16" (40 cm), 17" (42.5 cm), 18" (45 cm)				
Seat height (rear)	17" (43 cm)				
Seat height (front)	18.5" (47cm)				
Seat angle	4°				
Seat depth	15.75" (40 cm)				
Total length	36.6" (93 cm)				
Total height	25" - 33" (64-84 cm)				
Total width	22.5" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)	28" (71 cm)
Max User weight	220 lb (100 kg)			330 lb (150 kg)	
Total Weight	20.7 lb (9.4 kg)	20.92 lb (9.5 kg)	21.14 lb (9.6 kg)	21.8 lb (9.9 kg)	23.13 lb (5.5 kg)
Transfer Weight	10.13 lb (4.6 kg)	10.35 lb (4.7 kg)	10.57 lb (4.8 kg)	11.01 lb (5kg)	12.11 lb (5.5 kg)

PANTHERA S3 SWING SHORT SPECIFICATIONS

Seat width	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)
Backrest height	8" (20 cm), 10" (25 cm), 12" (30 cm), 14" (35 cm), 16" (40 cm), 17" (42.5 cm), 18" (45 cm)			
Seat height (rear)	17" (43 cm)			
Seat height (front)	18.5" (47 cm)			
Seat angle	4°			
Seat depth	13.75" (35 cm)			
Total length	30.7 - 35.4" (78-90 cm)			
Total height	25" - 33" (64-84 cm)			
Total width	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)	16.5" (42 cm)
Max User weight	220 lb (100 kg)			
Total Weight	20.1 lb (9.12 kg)	20.3 lb (9.21 kg)	20.49 lb (9.30 kg)	20.68 lb (9.39 kg)
Transfer Weight	9.42 lb (4.28 kg)	9.62 lb (4.36 kg)	9.81 lb (4.45 kg)	10.0 lb (4.54 kg)

panthera

Distributed exclusively
in North America by



Triumph Mobility Inc.

www.triumphmobility.com | info@triumphmobility.com | 1.855.546.0711
R_0716_S3SwingFamilyt-SalesSheet_EN | Printed in Canada