


# panthera<sup>®</sup>

RIGID WHEELCHAIRS



**TRIUMPH**  
MOBILITY

Exclusive Distributor in North America

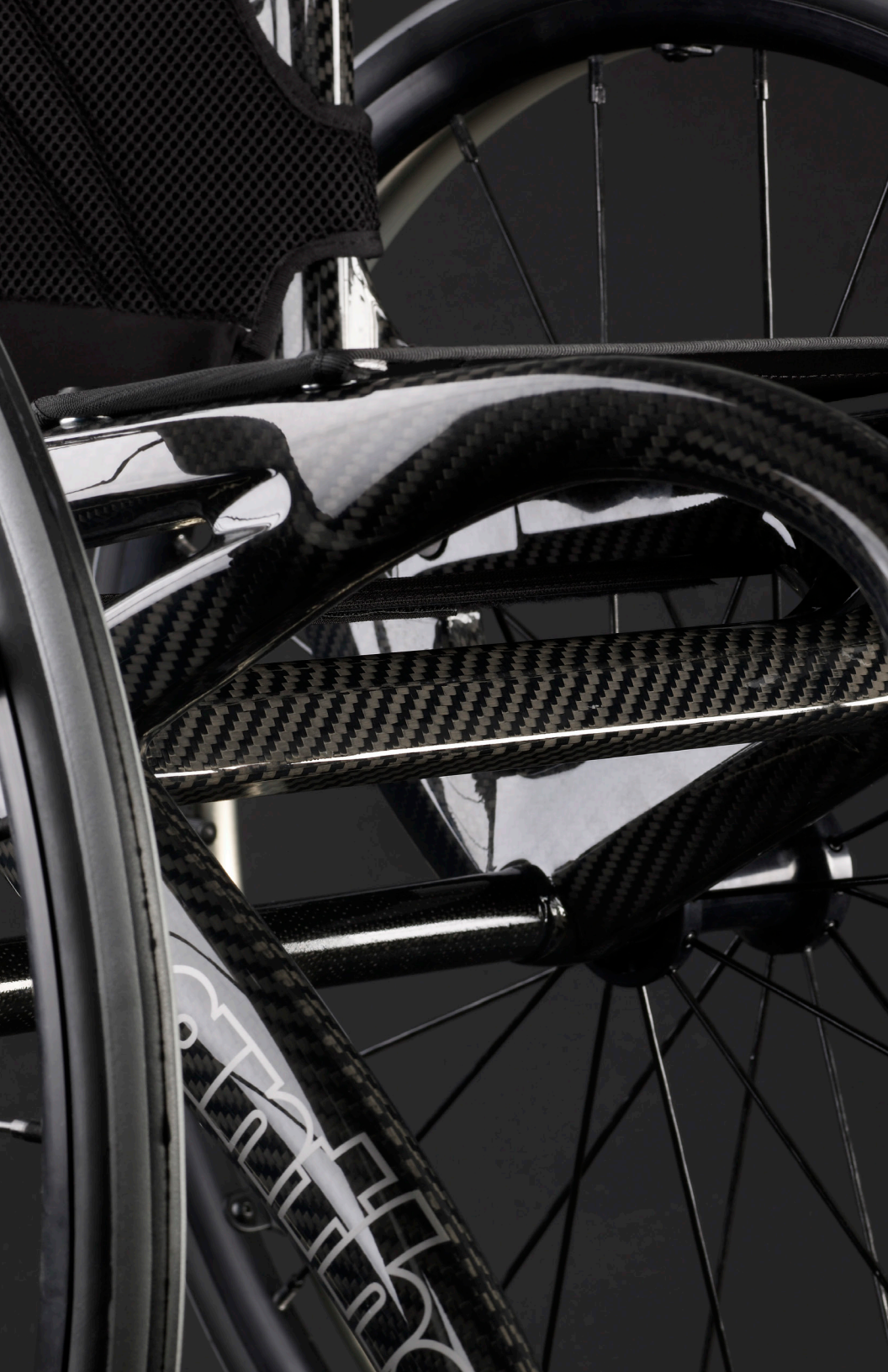

# panthera®

RIGID WHEELCHAIRS CATALOGUE

## CONTENTS

ABOUT TRIUMPH MOBILITY	5
PANTHERA RIGID WHEELCHAIRS	6
PANTHERA X	8
PANTHERA U3 LIGHT	10
PANTHERA U3	13
PANTHERA S3	16
PANTHERA S3 SWING	20
PANTHERA RIGID WHEELCHAIRS FOR CHILDREN	23
PANTHERA BAMBINO	24
PANTHERA BAMBINO HYBRID	26
PANTHERA MICRO	28



Triumph Mobility is driven to improve the daily lives of everyone we serve, helping them to move through their world with dignity and comfort.

Founded in 2012 with a goal to make and distribute innovative products that put the user first, Triumph is dedicated to helping clients live active and independent lifestyles. We are relentless in our efforts to source safe, dependable and ergonomic mobility solutions.

We provide mobility aids from rollators to wheelchairs, plus all of the accessories needed for these products. We also consult with clients to help them find the best products for their needs.

We manufacture the Escape Rollators and we are proud to be exclusive North American distributors for Panthera Rigid Wheelchairs from Sweden, and Rollz Motion Rollator/Transport Chairs from the Netherlands.

## VISION

Our vision is to become a global leader in providing dependable, comfortable, and innovative mobility and lifestyle solutions.

## MISSION

Triumph Mobility aims to improve the daily lives of everyone we serve, helping them to move through their world with dignity and comfort.

A Panthera makes most things in your daily life easier. It has excellent driving qualities, it is light to lift and it is comfortable to sit in.

This is because we have focused on the three most important features of a wheelchair: driving, sitting and lifting. These three things are done every day, all year round, and that is why they have to function in the best way possible.

## DRIVING

With a chassis of chrome molybdenum or carbon fibre and footbar and pushrims of titanium, the chair is extremely lightweight and incredibly easy to manoeuvre. A precisely calculated camber angle of 2.2° gives you great driving properties while the chair remains narrow.

## SITTING

A continuously adjustable angle for the backrest, luxuriously soft upholstery padded and adjustable at the top with a tightening strap throughout, allows you to form the back upholstery for exact support. The chairs have a 4° or 7° seat angle, and adjustable tightening strap in the seat.

This combination provides you with excellent options for finding a sitting position that is both comfortable and functional.

## LIFTING

Low weight, a low chassis and a natural, well-balanced grip make lifting the chair in and out of your car simplicity itself.



QUITE SIMPLY, THE BEST  
RIGID WHEELCHAIR.

panthera®

The right wheelchair gives you greater mobility, independence and can improve your health, setting you up for success to continue with your active lifestyle. Choosing the best model for you is an important decision and involves considering your needs, your mobility challenges and lifestyle factors.

When choosing a wheelchair, you will want to consider:

The activities that you pursue in your everyday life, whether you need to move over long distances, make frequent transfers in and out of the chair or to move the chair into storage in a car or elsewhere

- Your preference for manoeuvrability and your tolerance for how “tippy” the wheelchair might be (experienced users are more adept at handling a model without anti-tippers – the Panthera U3 Light model offers a lighter weight in exchange for a more tippy feel)
- Your preferences for footrest type
- The adjustability factors (many seats have different adjustability – the S3 models offer many options)
- Your preference for seating position—relaxed for comfort versus upright for a more active user (the S3 model offers a more relaxed seating position whereas the U3 model is more upright)
- Construction material (Ultralight vs. Light)
- Your size (short/large models and child versions such as the Bambino and Micro are available)

# panthera X™

BY FAR THE LIGHTEST AND MOST HIGH-TECH WHEELCHAIR EVER MADE,



Panthera X™ is by far the lightest and most high-tech wheelchair ever made, featuring the most advanced technology available today.

This elite wheelchair is designed for experienced active users who are able to handle a lightly balanced wheelchair without an anti-tip device and want to maintain an active lifestyle.

This wheelchair weighs in at an incredible 4.62 lb (without wheels) and features an X-shape side frame that makes it a dream to pick up and store in your car for transfer. Sleek mobility meets comfort with adjustability at both the backrest and footrest.

With chassis and backrest frame made of ultralight, ultra-strong carbon fibre material, the Panthera X offers an ultra-stable, dependable ride that relies on a rigid X frame for easy transfer in and out of your car. Panthera X has a built-in rear axle, and features 4 adjustments to its backrest at 1cm increments for added customization to change the centre of gravity. The backrest angle also has a continuously adjustable angle, and the footrest is adjustable for ultimate comfort.

## FEATURES

- Ultralight transfer weight\* frame starting at 4.62 lb (2.1 kg). With total weight: 9.68 lb (4.4 kg)
- Chassis and backrest frame both made from carbon fibre
- 24" Spinergy X rear wheel
- Single hand brake
- Rigid frame offers excellent stability
- Backrest with continuously adjustable angle and cushioned support for comfort
- 4 backrest positions provide high adjustability for frame balancing

\* 33 cm with X frame

## OPTIONAL ACCESSORIES

Panthera X™ offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Carbon Fibre Side Guards (each style sold separately)



Cushion X™



Transport Protection



Arm Hook



Tube Protection



Freewheel

## PANTHERA X SPECIFICATIONS

Seat width	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	17.7" (45 cm)
Backrest height	8.6 - 13.8" (22 - 35 cm)				
Seat height (rear)	17" (43 cm)				
Seat height (front)	18.5" (47 cm)				
Seat angle	7°				
Seat depth	13.8 - 18.1" (35 - 46 cm)				
Total length	32.5" (82.5 cm)				
Total height	25 - 33" (64 - 84 cm)				
Total width	21.25" (54 cm)	22.4" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)
Max user weight	220 lb (100 kg)				
Total Weight	9.7 lb (4.4 kg)	9.8 lb (4.45 kg)	9.9 lb (4.5 kg)	10 lb (4.55 kg)	10.1 lb (4.6 kg)
Transfer Weight	4.6 lb (2.1 kg)	4.7 lb (2.15 kg)	4.85 lb (2.2 kg)	4.96 lb (2.25 kg)	5.07 lb (2.3 kg)

# panthera® U3 LIGHT

AN ULTRA LIGHTWEIGHT RIGID WHEELCHAIR FOR EXPERIENCED USERS



The Panthera U3 Light is a high-end, ultra-lightweight wheelchair for the experienced active user who can handle a 'tippy' wheelchair without an anti-tip device.

The chair trades off some stability in exchange for lighter weight (transfer weight starts at 7.92 lb / 3.60 kg with wheels and sideguards removed) that makes driving the chair and transferring in and out of the car easier, without reducing the ergonomic features offered by the U3.


Two configurations offer further choice: the U3 Light with normal balance, and the "U3 Light L" with "tippy" balance (with the rear axle moved forward an additional 7/8" / 22 mm) and an adjustable backrest. The design is optimized for easy transfer in and out of vehicles, and features extremely good driving characteristics.

## FEATURES


- Lightweight frame: starting at 7.92 lb / 3.60 kg
- Tapered front and steeper footrest angle than the Panthera S3 model
- 7° seating angle and adjustable backrest
- Minimalist design
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- 24" Spinergy X rear wheels and 3" front casters designed for improved handling
- Stable, fixed frame for stability and adjustable balance point
- Single handed brake for agile movement

## CHASSIS COLOURS

The Panthera U3 Light comes standard in a pearl white or anodic black frame. It is also available in a series of premium colours which can be ordered at an extra cost.

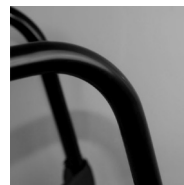

Anodic Black

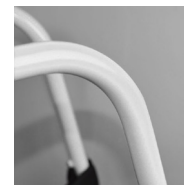


Pearl White




Classic Blue



Mean Green

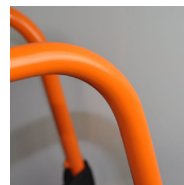

Panorange

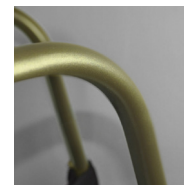

Brave Bronze

## OPTIONAL ACCESSORIES

Panthera U3 Light offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Height Adjustable Armrest



Carbon Fibre Side Guards



Foldable Push Handle



Seat Bag



Cushion



Freewheel

## U3 LIGHT SPECIFICATIONS

Seat width	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	17.7" (45 cm)
Backrest height	9.8" (25 cm), 11.8" (30 cm), 13.8" (35 cm)				
Seat height (rear)	17" (43 cm)				
Seat height (front)	18.5" (47 cm)				
Seat angle	7°				
Seat depth	13.8 - 18.1" (35 - 46 cm)				
Total length	32.5" (82.5 cm)				
Total height	25 - 33" (64 - 84 cm)				
Total width	21.25" (54 cm)	22.4" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)
Max user weight	220 lb (100 kg)				
Total Weight	13.2 lb (6.0 kg)	13.34 lb (6.05 kg)	13.45 lb (6.1 kg)	13.56 lb (6.15 kg)	13.7 lb (6.20 kg)
Transfer Weight	8.16 lb (3.70 kg)	8.3 lb (3.75 kg)	8.4 lb (3.80 kg)	8.5 lb (3.85 kg)	8.6 lb (3.90 kg)



# panthera U3

A HIGH-END RIGID WHEELCHAIR FOR THE ACTIVE USER



The Panthera U3 is a lightweight wheelchair designed for your active lifestyle. This chair has a tapered front and a steeper footrest angle that provides a more upright and active driving position compared to the Panthera S3 model. The Panthera U3 offers a stable, fixed frame and casters (wheels) with a specially designed tread make it very easy to manoeuvre. The 7° seating angle and adjustable backrest provide good seating comfort.

The chair's balance point can be adjusted to provide a feeling of security. Accessories such as armrests and side guards allow for even more customization, including adjustments to the height and positioning of the armrests that users can make themselves without tools.

You can also fold down the soft upper section of the side guards for easy movement in or out of the wheelchair. An attractively integrated anti-tip device can be folded in and out by the wheelchair user while seated.

## FEATURES

- Lightweight frame: starting at 9.85 lb (4.48 kg)
- Tapered front and steeper footrest angle than the Panthera S3 model
- 7° seating angle and adjustable backrest
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- Integrated anti-tip device can be folded in and out by the user
- 24" standard rear wheels and 3" front casters designed for improved handling
- Stable, fixed frame for stability and adjustable balance point
- Accessories allow for extensive customization and adjustments by the user




The Panthera U3 is recognized by the Canadian Association of Occupational Therapists.


# panthera U3

## CHASSIS COLOURS


The Panthera U3 comes standard in a beautiful anodic black frame. It is also available in a series of premium colours which can be ordered at an extra cost.

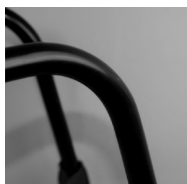

Anodic Black

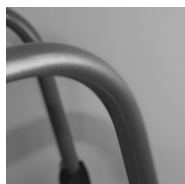


Grey Metallic

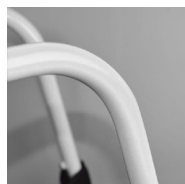


Pearl White




Classic Blue

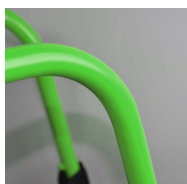

Mean Green

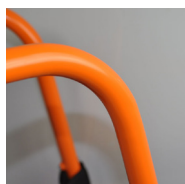

Panorange

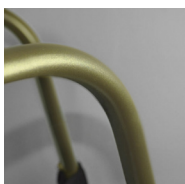

Brave Bronze

## OPTIONAL ACCESSORIES

Panthera U3 offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Spoke Protectors



Height Adjustable Push Handle



Foldable Push Handle



Push Bar



Height Adjustable Armrest



Carbon Fibre Side Guards



Anti-tip Device



Crutch Holder

## U3 SPECIFICATIONS

Seat width	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	17.7" (45 cm)
Backrest height	8" - 18" (20 - 45 cm)				
Seat height (rear)	17" (43 cm)				
Seat height (front)	18.5" (47 cm)				
Seat angle	7°				
Seat depth	13.75 - 18.1" (35 - 46 cm)				
Total length	33" (84 cm)				
Total height	25 - 33" (64 - 84 cm)				
Total width	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)
Max user weight	220lb (100 kg)				330 lb (150 kg)
Total Weight	17.7 lb (8.043 kg)	17.9 lb (8.13 kg)	18 lb (8.18 kg)	18.1 lb (8.22 kg)	18.8 lb (8.52 kg)
Transfer Weight	10 lb (4.48 kg)	10 lb (4.57 kg)	10.3 lb (4.66 kg)	10.5 lb (4.75 kg)	11 lb (4.97 kg)





# panthera s3

A LIGHTWEIGHT RIGID WHEELCHAIR FOR THE ACTIVE USERS



The Panthera S3 is a lightweight wheelchair designed for your active lifestyle. A stable, full-width frame and casters (wheels) with a specially designed tread make it very easy to manoeuvre. This chair provides a comfortable 7° seating angle with a right knee angle, a foot bar positioned further back and an adjustable backrest for a more upright, ergonomic sitting position that offers a feeling of security (see our Panthera U3 model if you are looking for a more upright, active driving position).

Accessories such as armrests and side guards allow for even more customization, including adjustments to the height and positioning of the armrests that users can make themselves without tools. You can also fold down the soft upper section of the side guards for ease transferring into or out of the wheelchair. An attractively integrated anti-tip device can be folded in and out by the user while seated.

## FEATURES


- Lightweight chassis: starting at 10.23 lb (4.65 kg)
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- 24" standard rear wheels and 5" front casters designed for improved handling
- Balancing continuously adjustable
- Backrest with continuously adjustable angle
- Adjustable backrest and seat upholstery
- Accessories allow for extensive customization and adjustments
- Integrated anti-tip device can be folded in and out by the user




The Panthera S3 Family of wheelchairs is recognized by the Canadian Association of Occupational Therapists.

## CHASSIS COLOURS

The Panthera S3 comes standard in a beautiful anodic black or grey metallic frame. It is also available in a series of premium colours which can be ordered at an extra cost.

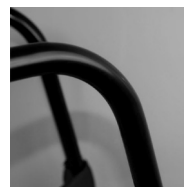

Anodic Black

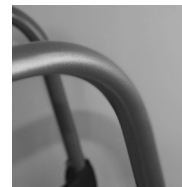


Grey Metallic




Pearl White




Classic Blue

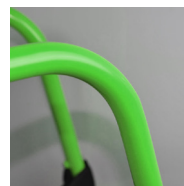

Mean Green

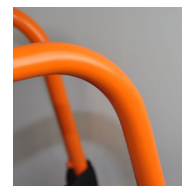

Panorange

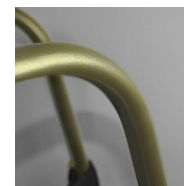

Brave Bronze

## OPTIONAL ACCESSORIES

Panthera S3 offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Spoke Protectors



Height Adjustable Push Handle



Foldable Push Handle



Push Bar




Height Adjustable Armrest



Carbon Fibre Side Guards



Anti-tip Device

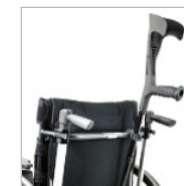

Crutch Holder

### S3 SPECIFICATIONS

<b>Seat width</b>	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	17.7" (45 cm)	20" (50 cm)
<b>Backrest height</b>	8 - 18" (20 - 45 cm)					
<b>Seat height (rear)</b>	17" (43 cm)					
<b>Seat height (front)</b>	18.5" (47 cm)					
<b>Seat angle</b>	7°					
<b>Seat depth</b>	15.75" (40 cm)					
<b>Total length</b>	30.7 - 35.5" (78 - 90 cm)					
<b>Total height</b>	25 - 33" (64 - 84 cm)					
<b>Total width</b>	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)	28" (71 cm)
<b>Max user weight</b>	220 lb (100 kg)			330 lb (150 kg)		
<b>Total weight</b>	18 lb (8.20 kg)	18.2 lb (8.28 kg)	18.4 lb (8.36 kg)	18.6 lb (8.44 kg)	19.16 lb (8.69 kg)	19.4 lb (8.78 kg)
<b>Transfer weight</b>	10.25 lb (4.65 kg)	10.43 lb (4.73 kg)	10.6 lb (4.81 kg)	10.8 lb (4.90 kg)	11.3 lb (5.13 kg)	11.5 lb (5.22 kg)

### S3 SHORT SPECIFICATIONS

<b>Seat width</b>	11.8" (30 cm)	13" (33 cm)	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	17.7" (45 cm)
<b>Backrest height</b>	8 - 18" (20 - 45 cm)					
<b>Seat height (rear)</b>	17" (43 cm)					
<b>Seat height (front)</b>	18.5" (47 cm)					
<b>Seat angle</b>	7°					
<b>Seat depth</b>	10.6 - 13" (27 - 33 cm)		13.8" (35 cm)			
<b>Total length</b>	28.75 - 33.5" (73 - 85 cm)					
<b>Total height</b>	25 - 33" (64 - 84 cm)					
<b>Total width</b>	20" (51 cm)	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)
<b>Max user weight</b>	220 lb (100 kg)			330 lb (150 kg)		
<b>Total weight</b>	17.7 lb (8.05 kg)	17.89 lb (8.13 kg)	18.08 lb (8.22 kg)	18.28 lb (8.31 kg)	18.48 lb (8.40 kg)	18.96 lb (8.62 kg)
<b>Transfer weight</b>	9.88 lb (4.49 kg)	10.05 lb (4.57 kg)	10.25 lb (4.66 kg)	10.45 lb (4.75 kg)	10.65 lb (4.84 kg)	10.9 lb (4.96 kg)

### S3 LARGE SPECIFICATIONS

<b>Seat width</b>	15" (39 cm)	16.5" (42 cm)	17.7" (45 cm)	20" (50 cm)
<b>Backrest height</b>	8" - 18" (20 - 45 cm)			
<b>Seat height (rear)</b>	18" (45.5 cm)			
<b>Seat height (front)</b>	19.5" (49.5 cm)			
<b>Seat angle</b>	7°			
<b>Seat depth</b>	17.7" (45 cm)			
<b>Total length</b>	30.7 - 35.4" (78 - 90 cm)			
<b>Total height</b>	26 - 34" (66.5 - 86.5 cm)			
<b>Total width</b>	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)	28" (71 cm)
<b>Max User weight</b>	330 lb (150 kg)			
<b>Total weight</b>	20.2 lb (9.17 kg)	20.3 lb (9.23 kg)	20.9 lb (9.47 kg)	21 lb (9.52 kg)
<b>Transfer weight</b>	12.4 lb (5.6 kg)	12.5 lb (5.67 kg)	13 lb (5.91 kg)	13.2 lb (5.98 kg)

### S3 SHORT LOW SPECIFICATIONS

<b>Seat width</b>	13" (33 cm)	14" (36 cm)	15" (39 cm)
<b>Backrest height</b>	8" - 18" (20 - 45 cm)		
<b>Seat height (rear)</b>	16" (40.5 cm)		
<b>Seat height (front)</b>	17.5" (44.5 cm)		
<b>Seat angle</b>	7°		
<b>Seat depth</b>	10.6 - 13" (27 - 33 cm)	13.8" (35 cm)	
<b>Total length</b>	28.75 - 33.5" (73 - 85 cm)		
<b>Total height</b>	25 - 33" (64 - 84 cm)		
<b>Total width</b>	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)
<b>Max User weight</b>	220 lb ( 100 kg)		
<b>Total weight</b>	17.9 lb (8.12 kg)	18.1 lb (8.21 kg)	19.3 lb (8.3 kg)
<b>Transfer weight</b>	10 lb (4.56 kg)	10.25 lb (4.65 kg)	10.45 lb (4.74 kg)

# panthera S3 SWING

A LIGHTWEIGHT RIGID WHEELCHAIR WITH SWING AWAY FOOTRESTS



The Panthera S3 Swing Wheelchair offers the same excellent features of the S3 Wheelchairs, but with divided swing-away footrests to allow for easy movement in and out of the chair.

The Panthera S3 Swing is a lightweight wheelchair designed for your active lifestyle. With a seat angle of 4°, it opens up new possibilities to find a good sitting position, even for those who use their legs to move about with the wheelchair, or for other reasons need to have a more flat seat angle.

The chassis seating angle, along with the adjustable backrest and seat upholstery and the backrest tilt function combine to provide extremely good seating comfort. A stable, fixed frame and casters (wheels) with a specially-designed tread make it very easy to manoeuvre. The chair's balance point can be adjusted to provide a feeling of security for every individual.

Accessories such as armrests and side guards allow for even more customization, including adjustments to the height and positioning of the armrests that users can make themselves without tools. You can also fold down the soft upper section of the side guards for easy movement in or out of the wheelchair.

The attractively integrated anti-tip device can be folded in and out by the wheelchair user while seated.

This model also has a divided swing-away footrest, with a user-friendly, discrete locking mechanism integrated inside the frame.

## FEATURES

- Lightweight frame: starting at 9.15 lb (4.160 kg)
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- 24" standard rear wheels and 5" front casters designed for improved handling
- Folding swing-away footrests
- 4° seat angle
- Integrated anti-tip device can be folded in and out by the user
- Balancing continuously adjustable
- Backrest with continuously adjustable angle
- Adjustable backrest and seat upholstery
- Accessories that allow for extensive customization and adjustments



The Panthera S3 Swing is recognized by the Canadian Association of Occupational Therapists.

## CHASSIS COLOURS

The Panthera S3 Swing comes standard in a beautiful anodic black frame. It is also available in a series of premium colours which can be ordered at an extra cost.



Anodic Black



Grey Metallic



Pearl White




Classic Blue



Mean Green



Panorange

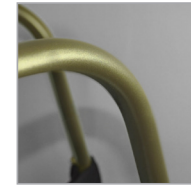

Brave Bronze

## OPTIONAL ACCESSORIES

Panthera S3 Swing offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Spoke Protectors



Height Adjustable Push Handle



Foldable Push Handle



Push Bar



Height Adjustable Armrest



Carbon Fibre Side Guards



Anti-tip Device



Crutch Holder

### S3 SWING SPECIFICATIONS

Seat width	14" (36 cm)	15" (39 cm)	16.5" (42 cm)	18" (45 cm)	20" (50 cm)
Backrest height	8 - 18" (20 - 45 cm)				
Seat height (rear)	17" (43 cm)				
Seat height (front)	18" (45 cm)				
Seat angle	4°				
Seat depth	15.75" (40 cm)				
Total length	36.6" (93 cm)				
Total height	25 - 33" (64 - 84 cm)				
Total width	22.5" (57 cm)	23.6" (60 cm)	24.8" (63 cm)	26" (66 cm)	28" (71 cm)
Max user weight	220 lb (100 kg)			330 lb (150 kg)	
Total weight	20.72 lb (9.40 kg)	20.9 lb (9.50 kg)	21.2 lb (9.60 kg)	21.8 lb (9.90 kg)	22.7 lb (10.30 kg)
Transfer weight	10 lb (4.60 kg)	10.4 lb (4.70 kg)	10.6 lb (4.80 kg)	11 lb (5.0 kg)	12 lb (5.5 kg)

### S3 SWING SHORT SPECIFICATIONS

Seat width	13" (33 cm)	14" (36 cm)	15" (39 cm)
Backrest height	8 - 18" (20 - 45 cm)		
Seat height (rear)	17" (43 cm)		
Seat height (front)	18" (45 cm)		
Seat angle	4°		
Seat depth	13.8" (35 cm)		
Total length	30.7 - 35.5" (78 - 90 cm)		
Total height	25 - 33" (64 - 84 cm)		
Total width	21.25" (54 cm)	22.5" (57 cm)	23.6" (60 cm)
Max User weight	220 lb (100 kg)		
Total weight	19.6 lb (8.9 kg)	19.8 lb (9.0 kg)	20 lb (9.1 kg)
Transfer weight	9.7 lb (4.4 kg)	9.9 lb (4.5 kg)	10 lb (4.6 kg)



QUITE SIMPLY,  
THE BEST RIGID  
WHEELCHAIRS  
FOR ACTIVE KIDS

**panthera**<sup>®</sup>

Finding the right wheelchair for a child can make a real difference in increasing their confidence and offering them the support they need to explore their world and keep up with their siblings and friends.

When choosing wheelchairs for a child, you'll want to consider your child's size, abilities and needs. Wheelchairs designed for children offer many adjustments so they can be adapted to their growth.

When choosing a child's wheelchair, consider:

- The environment that your child moves through and activities they pursue
- Adjustments available to meet their growth (many adjustments are available, including footrest height, angle and depth)
- Your child's size (the Bambino is a good everyday chair for children ages 4 to 8 and the Micro meets the needs of smaller children ages 1 to 5)

# panthera BAMBINO

A VERSATILE WHEELCHAIR FOR ACTIVE CHILDREN



The Panthera Bambino offers an easy-to-drive, everyday chair for children with many different needs. The chair offers your child a smooth, stable ride that can raise their confidence about their own mobility and their ability to engage with their busy world. The chair is ideal for children ages 4 to 12.

Like the adult versions, the Bambino provides a smooth ride, a wide range of accessories, and adjustment options that allow it to grow with the user. The chair is easy to drive, adapt, and lift, and an abduction front allows flexible leg positioning, offering 6 cm additional width at the knees with the same seat depth.

## FEATURES

- Lightweight frame: starting at 10.76 lb (4.89 kg)
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- Rear wheels from 20" to 24" and 5" front caster wheel
- Footrest adjustable in height, angle and depth
- Stable, fixed frame for stability and adjustable balance point
- Dual adjustable anti-tippers

## OPTIONAL ACCESSORIES

Panthera Bambino offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Spoke Protectors



Height Adjustable Push Handle



Push Bar



Side Guards



Cushion

## CHASSIS COLOURS

The Panthera Bambino comes standard in a beautiful anodic black or red frame. It is also available in a series of premium colours which can be ordered at an extra cost.



Anodic Black



Red



Grey Metallic



Classic Blue




Pearl White



Mean Green



Panorange

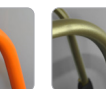

Brave Bronze

## BAMBINO SPECIFICATIONS

Seat width	9.45" (24 cm)	10.63" (27 cm)	11.8" (30 cm)	13" (33 cm)
Backrest height	6.7 - 9.8" (17 - 25 cm) or 10.6 - 13.8" (27 - 35 cm)			
Seat height (rear)	15" (39 cm)			
Seat height (front)	16.1" (41 cm)			
Seat angle	4°			
Seat depth	7 - 9.8" (18 - 25 cm)	9 - 11.8" (23 - 30 cm)		
Total length	25.6 - 28.1" (65 - 71.5 cm)			
Total height	22.6 - 25.6" (57.5 - 65 cm)			
Total width	20.5" (52 cm)	21.65" (55 cm)	22.8" (58 cm)	24" (61 cm)
Max User weight	132 lb (60 kg)			
Total weight	17.9 lb (8.12 kg)	18.4 lb (8.35 kg)	18.9 lb (8.58 kg)	19.4 lb (8.80 kg)
Transfer weight	10.8 lb (4.89 kg)	11.3 lb (5.14 kg)	11.8 lb (5.37 kg)	12.3 lb (5.6 kg)

# panthera BAMBINO HYBRID

A VERSATILE WHEELCHAIR FOR ACTIVE CHILDREN



The Panthera Bambino offers an easy-to-drive, everyday chair for children with many different needs. The chair offers your child a smooth, stable ride that can raise their confidence about their own mobility and their ability to engage with their busy world. The chair is ideal for children ages 4 to 12.

Like the adult versions, the Bambino provides a smooth ride, a wide range of accessories, and adjustment options that allow it to grow with the user. The chair is easy to drive, adapt, and lift, and an abduction front allows flexible leg positioning, offering 6 cm additional width at the knees with the same seat depth.

## FEATURES

- Lightweight frame: starting at 10.76 lb (4.89 kg)
- Chassis made from chrome molybdenum steel, rear axle from carbon fibre
- Rear wheels from 16" to 18" and 3" front caster wheel
- Footrest adjustable in height, angle and depth
- Stable, fixed frame for stability and adjustable balance point
- Dual adjustable anti-tippers

## OPTIONAL ACCESSORIES

Panthera Bambino Hybrid offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Spoke Protectors



Height Adjustable Push Handle



Push Bar



Side Guards



Cushion

## CHASSIS COLOURS

The Panthera Bambino Hybrid comes standard in a beautiful anodic black or red frame. It is also available in a series of premium colours which can be ordered at an extra cost.



Red



Anodic Black



Grey Metallic



Classic Blue



Pearl White



Mean Green



Panorange



Brave Bronze

## BAMBINO HYBRID SPECIFICATIONS

Seat width	9.45" (24 cm)	10.63" (27 cm)	11.8" (30 cm)
Backrest height	6.7 - 9.8" (17 - 25 cm) or 10.6 - 13.8" (27 - 35 cm)		
Seat height (rear)	13.9" (35 cm)		
Seat height (front)	15" (39 cm)		
Seat angle	4°		
Seat depth	7 - 9.8" (18 - 25 cm)	9 - 11.8" (23 - 30 cm)	
Total length	24.75" (63 cm)		
Total height	22.6 - 25.6" (57.5 - 65 cm)		
Total width	17.3" (44 cm)	18.5" (47 cm)	19.7" (50 cm)
Max User weight	132 lb (60 kg)		
Total weight	13.9 lb (6.31 kg)	14.4 lb (6.50 kg)	14.9 lb (6.77 kg)
Transfer weight	10.8 lb (4.89 kg)	11.3 lb (5.14 kg)	11.8 lb (5.37 kg)

# panthera MICRO

A LIGHTWEIGHT RIGID WHEELCHAIR VERY YOUNG CHILDREN



The Panthera Micro is a wheelchair for very young children, designed to meet the needs of those between 1 and 5 years of age (the youngest user we have met was 8 months old). The chair is intended for indoor use, and is extremely easy to drive for little people who want to be on the move!

Like the adult versions, the Micro provides a smooth ride and a wide range of accessories. Seat depth adjustment options and growth pad accessories allow the chair to grow with the user. The chair is easy to drive, adapt and lift (weight starts at 5.15 lb / 2.342 kg), and offers 16" or 18" rear wheels and a 3" single front caster wheel for easy rotation.

## FEATURES

- Lightweight frame: starting at 5.15 lb (2.342 kg)
- Chassis made from chrome molybdenum steel, rear axle from titanium (10° camber)
- 16" or 18" rear wheels and 3" front caster wheel
- Height-adjustable footplate
- No pushrims
- Anti-tipper

## OPTIONAL ACCESSORIES

Panthera Micro offers a great variety of high-quality accessories that help enhance the wheelchair comfort and performance. These can be purchased separately.



Spoke Protectors



Pushbar MICRO



Extended Backrest



Side Guards



Cushion

## CHASSIS COLOURS

The Panthera Micro comes standard in a beautiful Classic Blue. It is also available in a series of premium colours which can be ordered at an extra cost.



Classic Blue



Mean Green



Panorange



Brave Bronze

## MICRO SPECIFICATIONS

Seat width	9.4" (24 cm)	10.6" (27 cm)	9.4" (24 cm) Long	10.6" (27 cm) Long
Backrest height	8.7" (22 cm)			
Seat height (rear)	12.2" (31 cm)			
Seat height (front)	12.2" (31 cm)			
Seat angle	0°			
Seat depth	5.8 - 7.8" (15 - 20 cm)		7.8 - 9.8" (20 - 25 cm)	
Total length	24.4" (62 cm)		28.5" (72.5 cm)	
Total height	20.25" (51.5 cm)		21.25" (54 cm)	
Total width	18" (46 cm)	19.25" (49 cm)	18" (46 cm)	19.5" (49 cm)
Max User weight	66.14 lb (30 kg)			
Total weight	7.76 lb (3.52 kg)	7.86 lb (3.56 kg)	8.34 lb (3.78 kg)	8.44 lb (3.83 kg)
Transfer weight	5.16 lb (2.34 kg)	5.3 lb (2.38 kg)	5.45 lb (2.47 kg)	5.56 lb (2.52 kg)



PANTHERA

PANTHERA



# panthera®



151 Bentley Ave., Unit 4  
Nepean ON K2E 6T7  
Canada

t 1.855.546.0711  
e [info@triumphmobility.com](mailto:info@triumphmobility.com)  
w [triumphmobility.com](http://triumphmobility.com)



Panthera Wheelchairs Catalogue V. 200615

The information on this document is accurate as of July 15, 2020. Triumph Mobility and Panthera Wheelchairs reserve the right to update any information without prior notice. For the most up-to-date version of this catalogue, please visit our website or contact Triumph Mobility.